



Oviyam 2.0 Setup

1. Copy Oviyam2.war file into your dcm4chee installation's deploy folder.
2. Point your browser URL to <http://your-dcm4chee-ip:8080/Oviyam2>
3. Use dcm4chee's login credentials for login.

4. **Settings:** The gear icon in the top left provides access to the setting options.

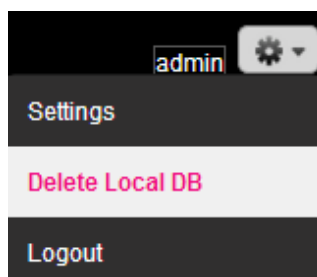

5. Setup dcm4chee's AE Configuration (Multiple WADO supported DICOM Servers can be configured). Only users with “WebAdmin” rights or higher can access this page.

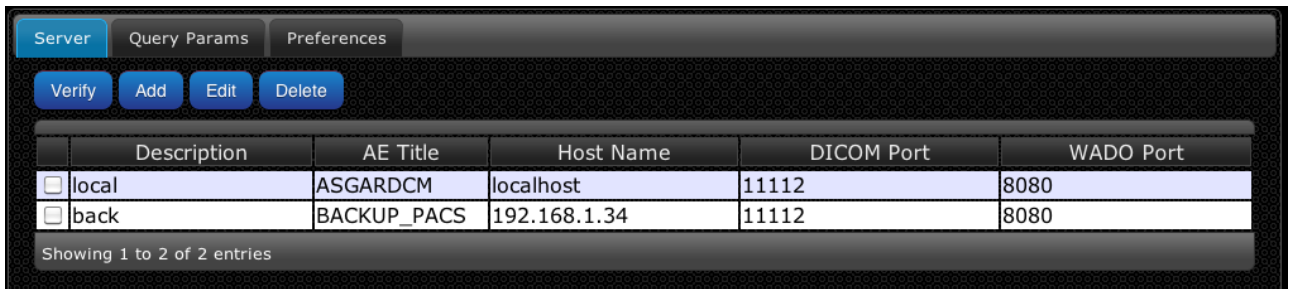


The screenshot shows a web interface for configuring DICOM servers. At the top, there are tabs for 'Server', 'Query Params', and 'Preferences'. Below the tabs are buttons for 'Verify', 'Add', 'Edit', and 'Delete'. The main area contains a table with the following data:

	Description	AE Title	Host Name	DICOM Port	WADO Port
<input type="checkbox"/>	local	ASGARDCM	localhost	11112	8080
<input type="checkbox"/>	back	BACKUP_PACS	192.168.1.34	11112	8080

Below the table, it says 'Showing 1 to 2 of 2 entries'.

6. Quick Search: Predefined Search filters can be created and assigned a name that appears as a button on the query screen. Examples:
- Today CT : (Query date/Current Date will be used as “Study Date” filter & “CT” will be used as “Modalities in Study” filter)
 - All CR : (“Modalities in Study” filter - “CR”)

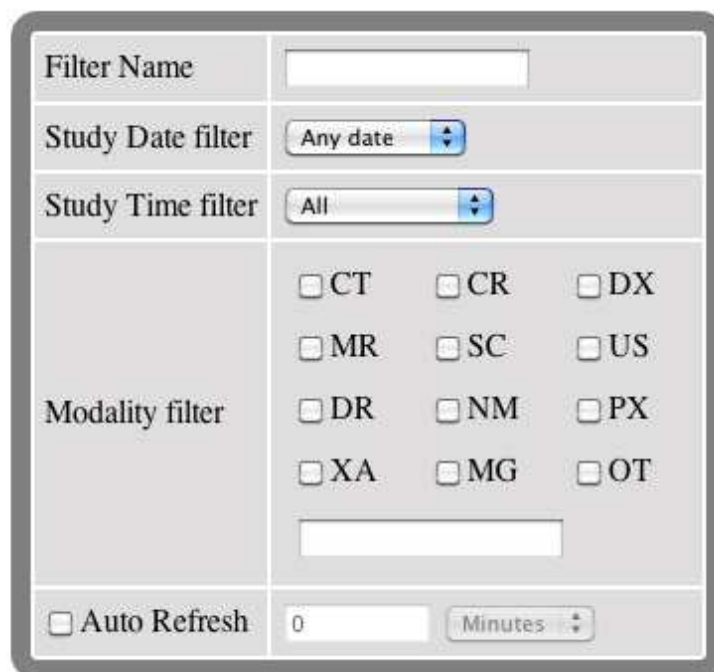

The screenshot shows a dialog box for configuring search filters. It has the following sections:

- Filter Name:** A text input field.
- Study Date filter:** A dropdown menu with 'Any date' selected.
- Study Time filter:** A dropdown menu with 'All' selected.
- Modality filter:** A grid of checkboxes for different modalities: CT, CR, DX, MR, SC, US, DR, NM, PX, XA, MG, OT. Below the grid is a text input field.
- Auto Refresh:** A checkbox and a numeric input field set to '0', followed by a 'Minutes' label and a dropdown arrow.

Quick Search filters will be shown as a buttons on the query screen.



7. Preferences

- Slider Visibility – User can opt whether the image slider should be shown or hidden in the image viewer. Default value is “Hide”.
- Session timeout – Default value is “30” mins.
- Theme – User can opt a UI theme from the list of themes. Default is “Dark Hive”.

8. Delete Local DB

Users can delete the studies downloaded in the local DB (browser). Local studies will be deleted, if the browser cache is cleared.

